

Invacare® | **Essential Plus**



Warranty

4 years foam
2 years cover



Cost effective pressure reducing mattress

The *Invacare Essential Plus* is an entry level pressure ulcer prevention mattress designed for the community healthcare environment. Practical and durable, this mattress offers comfort, support and pressure reduction at a cost effective price. The *Essential Plus* is suitable for those considered at 'High Risk' of developing pressure ulcers.



Yes, you can.®



Features and Options

Effective Pressure Ulcer Prevention

The core of the *Invacare Essential Plus* is constructed from a series of 'keyhole cut' cells, which move independently to evenly distribute weight, promoting greater comfort and facilitating pressure ulcer prevention.



Turn and rotate

Double sided castellated foam allows the mattress to be turned or rotated for increased product longevity.



Available in the following sizes:

L2000 x W855 x D140mm
Part no: 1551553

L1950 x W880 x D140mm
Part no: 1551551

Cover is resistant to fluid ingress


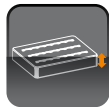



The outer cover is a multi stretch, water resistant, polyurethane fabric providing excellent recovery and durability which help to reduce shear and friction forces. The cover has a white inner liner allowing easy inspection for strikethrough or other cover damage.



Zip around three sides


Allows easy inspection of the foam core and easy removal for laundering.

Technical data

				
<i>Invacare Essential Plus</i> Mattress sizes : 2000 x 855 mm and 1950 x 880 mm	Height: 140 mm	Weight: 14 kg	Max. User Weight: 140 kg	Cleaning: 80°

Mattress is manufactured to:
■ BS7177:1996 Crib 5.

Colour

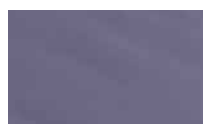

Baltic Blue



Invacare Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2012 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Essential Plus - UK - 12/2012 - 1556578



Yes, you can.®

1/ Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.