

Invacare® | *Etude® Medley Range*



3 year warranty



Cost effective patient care

The **Etude Medley** from Invacare offers the perfect balance of budget, functionality and quality requirements. The four-sectioned profiling mattress support allows for a range of height combinations, providing an ideal position for the carer and comfort for the user.

The latest addition to the **Etude Medley** range is the **Etude Medley Low** which lowers to a 22cm height to maximise safety for vulnerable clients who are prone to unexpected bed leaving. A specially designed bed end fascia then enables both the **Etude Medley** and **Etude Medley Low** to be transformed into the **Invacare Etude Medley Select**. This stylish bed is available in a beech finish that gives the type of homely feel that will integrate into most community and Residential Care environments. All beds are easy to transport and assemble.



Yes, you can.®

Improved functionality and efficiency

The *Invacare Etude Medley* has been specifically designed for quick transportation, to meet the needs of the transient working environments of Community Care provision. The **Medley** can be supplied on a transport kit so that the bed can easily be moved and installed. In addition, any assembling can be done without the use of tools or technicians.

Height is no longer an issue

The flexible design of the **Etude Medley** allows the bed to be mounted at a low height or a higher position, creating a range of working heights for the carer and transfer heights for the patient. This is facilitated by the unique feature of the bed end, which enables the mattress support to be mounted at both upper and lower positions.



Options and features



Double height positions on bed ends

Adjustable height range of 33 - 73 cm or 40 - 80 cm are available for the **Etude Medley**.



Four section profiling mattress support

Etude Medley comes with a full electric mattress support platform to suit individual patient requirements.



Full length side rails

The **Etude Medley** is available with wooden or steel full length side rails. Three quarter length metal side rails are also available in a range of heights.

Technical data


Etude Medley Width outside: 102 cm Width inside: 90 cm	Length outside: 222 cm Length inside: 208 cm (+ 20 cm extensions)	Height adjustable: 33 - 73 or 40 - 80 cm	Mattress support dimensions 4 sectioned: 75 - 35 - 34 - 56 cm	0 - 70° Backrest angle	Thigh angle: 0 - 29° Legrest lift: 0 - 20°

Etude Medley 11° Tilt option*	69 kg	SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

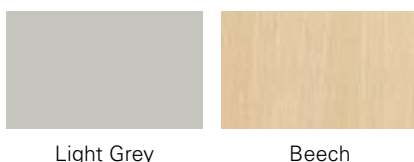

Light Grey

Beech

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk



Options and features

Low height to maximise safety

The height of the **Etude Medley Low** bed offers increased security for clients enabling them to sleep with real peace of mind. The **Etude Medley Low** sits to the floor, (22 cm in height) which combined with a floor mattress, means that if a fall does take place any harm is reduced. Carers also have the option to raise the bed to a more suitable height, when necessary.

The **Etude Medley Low** makes it easier for clients to get in and out of bed, providing more independence.

The **Medley Low** is dual height adjustable, offering a 22 - 62 and 29 - 69 cm height option.



Easy to transport

The **Etude Medley** and **Etude Medley Low** are easy to assemble and dismantle and can be transported using transporter kits.



Floor mattress helps to prevent injuries

In the event of a client falling from the bed, a foldable mattress can be placed on the floor and attached to the **Medley Low** for added safety.

The part number for the foldable mattress is 1494608.

Technical data

Etude Medley Low



Width outside: 102 cm
Width inside: 90 cm



Length outside: 222 cm
Length inside: 208 cm
(+ 20 cm extensions)



Height adjustable: 22 - 62 or 29 - 69 cm



Mattress support dimensions 4 sectioned: 75 - 35 - 34 - 56 cm



0 - 70° Backrest angle



Thigh angle: 0 - 29°
Legrest lift: 0 - 20°

Etude Medley Low



11° Tilt option*



69 kg



SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work



Light Grey



Beech

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Homely design

The **Etude Medley Select** offers a more attractive, homely design for **Etude Medley** and **Etude Medley Low**. Also available as a retro-fit, this aesthetic bed end cover fits with precision onto the bed actuators and can be removed at any time.

The leg supports are equipped with a Rastofix to enable individual positioning for the legs. The **Etude Medley Select** can be dismantled and stored easily using a transport adapter for moving and storage.




Options and features



Select design for Medley Low
The homely design of **Medley Select** is also available with the **Medley Low** model.



Lifting pole
Can be adjusted in both height and depth.

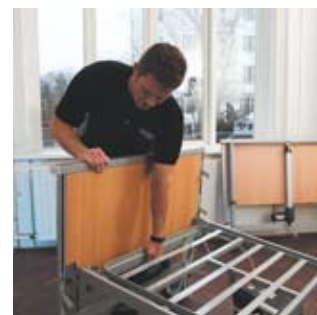

Medley can easily be dismantled by one person without using tools.

Technical data

Etude Medley Select	Width outside: 105 cm Width inside: 90 cm	Length outside: 225 cm Length inside: 208 cm (+ 20 cm extensions)	Height adjustable: 33 - 73 or 40 - 80 cm Medley Low: 22 - 62 or 29 - 69 cm	Mattress support dimensions 4 sectioned: 75 - 35 - 34 - 56 cm	0 - 70° Backrest angle	Thigh angle: 0 - 29° Legrest lift: 0 - 20°

Etude Medley Select	11° Tilt option*	76 kg SWL 180 kg (28.5 stone)

* Please notice that the tilt function is only practised for users under medical supervision.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

For mattress compatibility, please refer to the user manual which can be found at www.invacare.co.uk

Frame colour and wood work



Light Grey

Beech

INVACARE® Limited
Pencoed Technology Park - Pencoed - Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.

Etude Medley - UK - 09/2010
1484970



Yes, you can.®

