



*Invacare*<sup>®</sup> | **Matrix<sup>®</sup> Flo-tech<sup>™</sup> Contour Visco**

Contoured Visco memory foam that conforms to the shape of the body

**matrix**<sup>®</sup>  
SEATING SERIES

Trusted Solutions,  
Passionate People



The *Invacare* **Matrx Flo-tech Contour Visco** cushion is moulded from Visco-elastic foam. Its contoured surface area conforms to the body, improving posture, increasing comfort and maximising the surface area in contact with the user, thereby increasing pressure redistribution.

## Features and Options



### Deep moulded low memory foam cushion

The visco-elastic foam conforms to the body. By maximising the contact area it increases pressure redistribution and delivers exceptional levels of comfort. The deep foam makes this cushion particularly suitable for heavier and taller people.

### Ingress Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base with carry handle for easy transportation.

### Available in 51 sizes

The *Invacare* **Matrx Flo-tech Contour Visco** is available in a range of 51 sizes, from 36 x 33 to 51 x 46 cm (360 x 330 to 510 x 460 mm). (See size chart for details).

### Four foam densities available

Available in four densities to accommodate a range of client weights - Soft (up to 51 kg), Medium (up to 76 kg), Firm (up to 108 kg) and Very Firm (up to 127 kg).

### Contoured seating surface with pommel and side support

Offers increased lateral support and pelvic stability to enhance posture.

## Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see the [Invacare website](#).


	Various sizes	Height	Total product weight	Max. user weight <sup>1</sup>	Max. cleaning temperature	Warranty
<b>Matrx Flo-tech Contour Visco</b>	A range of 51 sizes	120 mm	1.6 kg	soft: 51kg (8 stone) medium: 76kg (12 stone) firm: 108kg (17 stone) very firm: 127kg (20 stone)	80°	3 years

**Disclaimer** 1. Guide maximum user weight limits are set assuming that the cushion selected is of appropriate width and depth for the intended user/s.


Manufactured to European Standards:  
Interior Foam: BS EN 1021-1, BS EN 1021-2, BS 5852: 1982 Crib 5

## Colours

Black Smooth

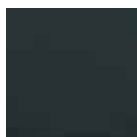

Black wipedown

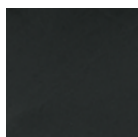

Black spacer fabric



Baltic Blue



Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Limited  
Pencoed Technology Park  
Pencoed - Bridgend CF35  
5AQ - United Kingdom

Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [ireland@invacare.com](mailto:ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

Invacare Ireland Ltd  
Unit 5 - Seatown  
Business Campus  
Seatown Road  
SWORDS - County  
Dublin - Ireland



**Yes, you can.®**

