

Invacare® | **I-Transia™** **Ceiling Hoist Solutions**

NEW



(250kg hoist option)



Invacare I-Transia - Solutions for Community and Residential Care

The *Invacare I-Transia* is a reliable and highly cost effective ceiling hoist that is perfect for community and residential care. Designed to meet the wide range of day-to-day lifting needs encountered by carers, the *I-Transia* allows transfers to be carried out quickly, safely and comfortably. Providing a range of simple solutions, such as; single track, room covering (XY) system and a free standing gantry, the *I-Transia* offers a solution for the most commonly found transfer issues.



Yes, you can.®



No nonsense installations

To ensure a fast and efficient service, **I-Transia** overhead hoisting systems have been designed to cover the most commonly requested formats.

1. Assess. 2 Order. 3. Install

Providing a simplified ordering process the **I-Transia**'s straightforward kit structure, significantly reduces the time from assessment to installation.

Peace of mind for the lifetime of the installation

The **I-Transia** ceiling hoist and track system is manufactured with an emphasis on technology and high quality specification. Supported by a network

of specially selected and highly skilled installers, you can be assured of safety and reliability throughout the life of the installation.

Features and Options



Easy charging – place the hand control in its holder




Snaphook for easy attachment of spreader bar



Emergency lowering and stop facility

Track Options

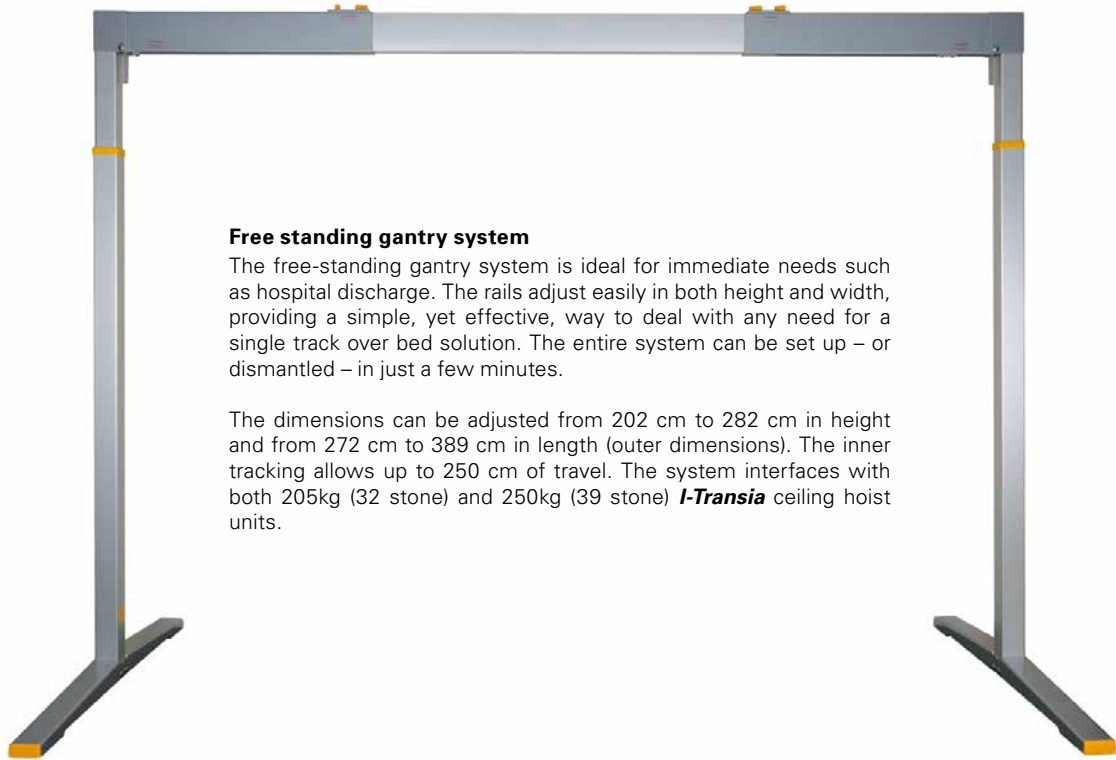

Free standing gantry system

The free-standing gantry system is ideal for immediate needs such as hospital discharge. The rails adjust easily in both height and width, providing a simple, yet effective, way to deal with any need for a single track over bed solution. The entire system can be set up – or dismantled – in just a few minutes.

The dimensions can be adjusted from 202 cm to 282 cm in height and from 272 cm to 389 cm in length (outer dimensions). The inner tracking allows up to 250 cm of travel. The system interfaces with both 205kg (32 stone) and 250kg (39 stone) **I-Transia** ceiling hoist units.



Single track system

The single-track system comprises a single rail which is fixed to the ceiling or a wall.

A single-track can be installed, up to 3 or 4 metres and takes the 205kg **I-Transia** ceiling hoists. 250kg hoists can be mounted on track up to 4 metres ceiling mounted and 3.5 metres if mounted to a wall.






Room covering (XY system) track system

This system allows transfers to be carried out in any area of the room. Ideal for when flexibility is required.


Technical data

						
<i>Invacare I-Transia</i> 205 kg	350 mm	194 mm	156 mm	2500 mm	40mm/sec at 85 kg 80mm/sec at max 5 kg	631 mm
<i>Invacare I-Transia</i> 250 kg	350 mm	194 mm	156 mm	2500 mm	40mm/sec at 85 kg 80mm/sec at max 5 kg	631 mm

		
<i>Invacare I-Transia</i> 205 kg	8 kg	205 kg (32 stone)
<i>Invacare I-Transia</i> 250 kg	8 kg	250 kg (39 stone)

Bariatric use

Colours

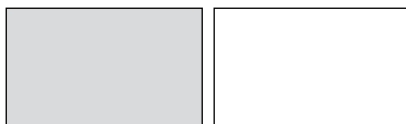

Rail colour: Grey Hoist colour: White

Lifts per charge; Aprox 55. (user weight 85 kg, 1000 mm up/down)

Invacare Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2013 Invacare International Sàrl
All rights reserved.

I-Transia - UK - 04/2013 - 1560319