

Invacare® | *ScanBed 750*



Excellent ergonomic features and comfort

The *Invacare ScanBed 750* excels in its ergonomic design. Not only does the backrest allow an anatomically correct return during elevation, the lower leg section can also be positioned beyond the side structure and seat section. In practical this means a smooth, effortless transfer from a bed to armchair position providing excellent comfort in the heart-/Fowler positioning.



Yes, you can.®



Freedom, independence & maximum comfort

Invacare **ScanBed 750** affords the user an exceptionally high degree of freedom and autonomy. Reclining back rest and adjustable thigh section facilitate maximum comfort in sitting & resting positions, and help prevent the occurrence of pressure sores.

A majestic lifting range of 35-82 cm makes getting in and out of bed trouble-free and supports easy transfer, for instance, from bed to wheelchair. User safety is further reinforced by the support handle.

Cleaning & maintenance - no problem!

Cleaning the **ScanBed 750** couldn't be easier! Both the backrest & leg rest fold away vertically affording easy access to the rest of the mattress support platform. Moreover, all mattress support electric components can be replaced without the user having to leave the bed.

Safety


Enjoy life with the **ScanBed 750**. Guaranteed for long-term use, the **ScanBed 750** offers excellent stability and good resistance against wear and tear. The results speak for themselves.

Designed to comply with all European Standards for Nursing Beds (EN 1970), the **ScanBed 750** passes with flying colours. During independent testing not a single adverse comment was made.


Features and Options




Invacare **Scala Basic 2**
Collapsible steel side rail. Compliant with the new standard for beds IEC 60601-2-52.

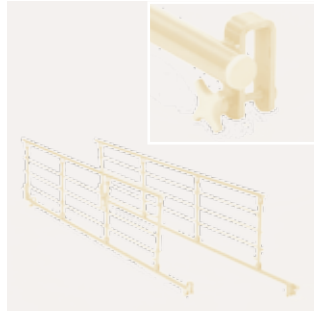

Invacare **Scala Basic Plus 2**
Collapsible steel side rail with plastic inserts at the mounting brackets. Compliant with the new standard for beds IEC 60601-2-52.

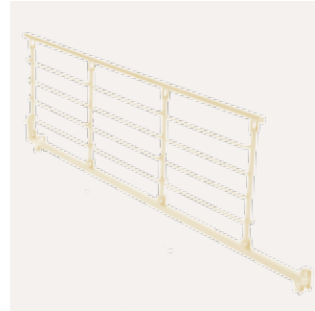


Invacare **Scala Medium 2**
Collapsible steel side rail for higher mattresses. Compliant with the new standard for beds IEC 60601-2-52.

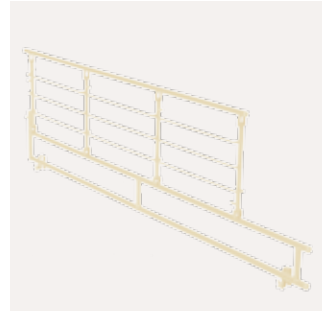

Invacare **Scala Decubi 2**
Collapsible steel side rail for extra high mattresses. Compliant with the new standard for beds IEC 60601-2-52.



Features and Options

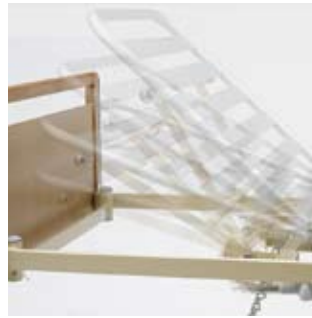

Auto regression
Increased comfort and reduced shear and friction when elevating the back rest.



Bed extension
Extension of the mattress support integrated.



Control Box
Cables and control box is coloured for easy mounting.



Snaplock function
Snaplock function separates mattress support from base.



Detachable bedends
Detachable bedends in a range of designs in beech or cherry.



Sanne bed ends
Available in beech, cherry, mix (cherry/beech) and in a high/low version.



Siderails
Siderails for wooden covering in standard and extended length in beech or cherry wood.




Hand control
The colours makes it easy to benefit from all the positioning options.



Support handles
Available in 5 different sizes make getting in and out of bed much easier for the patient. Is easy to place on the bed frame.



Central breaks
Available with a wheel diameter in 100, 125 or 150 mm.

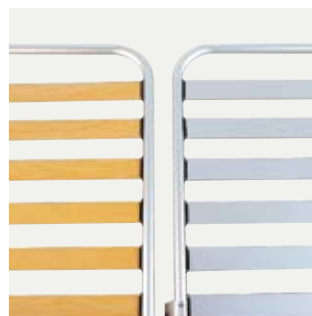

Mattress support
Mattress support with choice of aluminium or wooden slats.



Extra width
Beds with extra width in 105 cm or 120 cm comply with special needs. Mattress support in steel.

Please see *Invacare Bed Accessories* brochure for more accessories.



Technical data


SB 750	Width inside: 85 / 90 cm Width outside: 95 / 101 cm	Length outside: 213 cm Length inside: 208 cm Mattress: 200 cm	Height adjustable: 35-82 cm Total height: 160-207 cm	81 - 25 - 34 - 60 cm (4 - sectioned)	Under bed clearance: 17 cm	0° - 70° Backrest angle

SB 750	0° - 30° Thigh angle	0° - 17° Minimum angle below horizontal	0° - 16° Full legrest lift	32° (Angle below frame) 18° (Angle above frame)	SWL 220 kg Max. user weight 185 kg	Bariatric use*

*: Available in Heavy-duty version. Safe working load of 250 kg.

NB/ The bed must not be used by patients under 12 years of age, or by patients with body size equivalent to an average 12 year old or smaller.

Frame colour and wood work

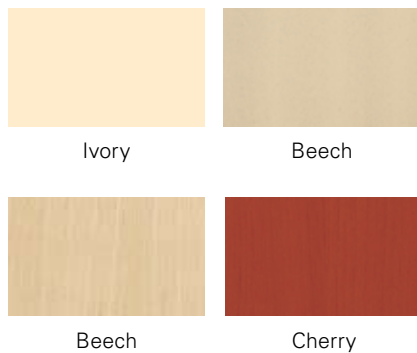

Invacare International Sarl
Route de Cité-Ouest 2 - 1196 Gland
Switzerland
Tel: +41 (0)22 354 60 10
Fax: +41 (0)22 354 60 11
hqeurope@invacare.com
www.invacare.eu.com

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.

SB 750 - EU - 08/2010



Yes, you can.®

