

Invacare® | **Softform Excel**



Warranty
Foam 2 Years
Cover 2 Years



Pressure reducing mattress

The *Invacare Softform Excel* provides a high level of comfort. The pressure relieving mattress core is constructed from high specification foam which provides even weight distribution for greater comfort and pressure re-distribution. With a robust foam base it prevents the mattress from bottoming out and is suitable for those considered 'At Risk' of developing pressure ulcers.



Yes, you can.®




Features and Options

Effective pressure ulcer prevention

The core of the *Invacare Softform Excel* is constructed from a series of 'keyhole cut' cells, which move independently to evenly distribute weight, promoting greater comfort and facilitating pressure ulcer prevention. The robust base features a firm layer of high density foam to provide substantial support.

Multi stretch PU waterproof cover

The cover fabric provides excellent recovery and durability, and helps to reduce shear and friction forces. The base of the mattress is a tough coated polyurethane (PU) which protects the mattress when used with electric profiling bed frames.

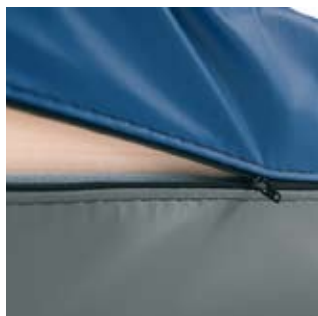

Zip around three sides

Allows easy inspection of the foam core.



Rotate-only mattress

Rotate-only turning cycle will increase product longevity. Non turn reduces the amount of manual handling required.



Strong, durable construction

The *Invacare Softform Excel* is an established choice of healthcare professionals.

Technical data

Softform Excel	Standard width: 80, 85.5, 88, 90 cm	Standard length: 195 cm	Height: 15.2 cm	Weight: 10 kg	Max. Client Weight: 140 kg (22 stone)

	Softform Excel	Cleaning: 80°
--	-----------------------	---------------

Manufactured to European Standards:
BS7177: Crib 5

INVACARE® Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be correct at time of print.
Invacare® reserves the right to alter product specifications without prior consultation.

Softform Excel - UK - 10/2010
1506260

Colour



Blue (PU)



Yes, you can.®