

Invacare® | *Softform Premier MaxiGlide*



Warranty
Foam 8 Years
Cover 4 Years

Designed to reduce shear and friction

The *Invacare Softform Premier MaxiGlide* mattress offers a new and innovative glide mechanism that significantly reduces shear and friction forces when used with an electric profiling bed. The patented design of the mattress allows the top surface to provide excellent pressure reducing properties along with a high level of comfort and support, whilst the base of the mattress uniquely separates to conform to the bedframe. The *Softform Premier MaxiGlide* is suitable for use with those who are considered to be at 'Very High Risk' of developing pressure ulcers.



Yes, you can.®



Effective pressure ulcer prevention

The **Sofform Premier MaxiGlide** features a deep-contoured insert pad constructed from 'keyhole' cut castellated foam cells, which move independently to maximise the body contact area and optimise pressure reduction.

Reduces the need for manual handling

The profiling mattress maximises the benefits of an electric profiling bed, delivering improved positioning whilst reducing the requirement for repositioning of the client.

Multi stretch cover

The mattress cover is multi stretch and conforms to the deep contoured castellations on the upper surface of the mattress, minimising rucking and further reducing the potential for tissue damage.

Features



Smooth sliding surface



Top surface remains with patient




Bottom surface profiles with the bedframe

Features



Patented design

The "Glide" mattress system is only available from Invacare and approved distributors.

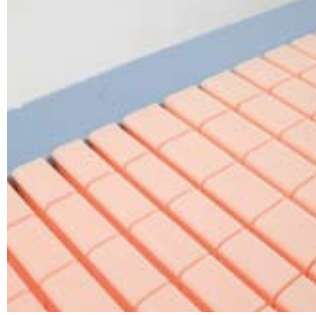

Modular Construction

Enables easy replacement of all parts, increasing product longevity.



Pressure Ulcer Prevention

The mattress provides an excellent pressure reducing surface, designed to maximise the body contact area and prevent the build up of pressure ulcers*.



Cover is resistant to fluid ingress

Features include two-way stretch, vapour-permeable and waterproof polyurethane cover with high frequency welded seams and a full length flap concealing the zip.



Fully detachable top cover and toughened polyurethane base

The toughened PU coated base protects the mattress against metal bed frames. The top cover is easily removed for laundering and foam inspection.



Does not require turning or rotating


Mattress does not require turning or rotating, significantly reducing manual handling risk.

Technical data


**Sofform
Premier
MaxiGlide**



Std Width:
80, 85.5, 88, 90 cm

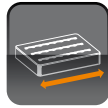

Std Length:
195-210 cm

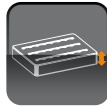

Height:
15.2 cm



Weight:
13.5 kg



Max. Client Weight:
247.6kg / 39 stone



Cleaning: 80° C
Refer to
instructions
for use and relevant
infection control
guidelines

Manufactured to European Standards:
Interior Foam: BS 5852 Crib 5
Cover: BS 7175 Crib 5

Colour



Baltic Blue

INVACARE® Limited
Pencoed Technology Park - Pencoed -
Bridgend CF35 5HZ
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
www.invacare.co.uk
Sales Order E-mail:
ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS - County
Dublin
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2010 Invacare® International Sarl
All rights reserved.

All information quoted is believed to be
correct at time of print.
Invacare® reserves the right to alter
product specifications without prior
consultation.

Sofform Premier MaxiGlide - UK
11/2010 - 1538670

* Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. A range of assessment scales can be used as a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.



Yes, you can.®