

# Invacare® | Storm<sup>4</sup> Max



## Much more than just a larger power chair!

The new *Invacare Storm<sup>4</sup> Max* is designed specifically for the needs of clients with a larger bodyshape. Its longer chassis supports the shifted centre of bodymass for improved weight distribution and swift mobility. Extra strong arm- and legrest support the clients needs to push themselves back to re-position or do some type of pressure relief.

The **G-Trac<sup>®</sup>** option is also available, realizing the ultimate driving performance and traction control.

The stylish **Storm<sup>4</sup> Max** completes the **Storm<sup>4</sup>** family and conserves the trouble-free servicing.



**Yes, you can.®**



**Features and Options**

**Tilt and recline**

For individual comfort and positioning the **Storm<sup>4</sup> Max** offers powered tilt and recline options. The backrest recliner has a range of 0° - 25° with a high pivot point in order to reduce shearing forces.

**Ease of use**

Just as the entire *Invacare Storm<sup>4</sup>* family is characterised by its ease of use, so is the **Max**. Batteries and electronics are easily accessible, making servicing simple. The depth of the seat can be simply adjusted in 3 steps which are indicated on the seat frame. An adapter to further reduce the seat depth is offered as an option.



**Armrest**

The strong and stylish armrest can be mounted in several positions as well as adjusted in width.



**Footboard**

Flip-up footboard supports various positions of the feet. It's strong enough to allow the clients to push themselves back to reposition or do some type of pressure relief.




**Calf pads**

The extra wide Max calf pads are individually adjustable in depth and angle to give optimum support.

**Technical data**

<b>Storm<sup>4</sup> Max</b>	490 - 590 mm 550 - 650 mm 610 - 710 mm	560, 595, 630 mm <sup>1</sup>	480, 510 mm	570 mm <sup>2,3</sup> , 630 mm <sup>3</sup>	360 - 480 mm	0°- 16° <sup>4</sup> , 0°- 9° <sup>5</sup>	0 - 25°	1390 mm
<b>Storm<sup>4</sup> Max</b>	2080 mm	60 mm	6°(10,5 %)	350W	31 km <sup>6</sup>	6 km/h	220 kg	

**Standard colours**

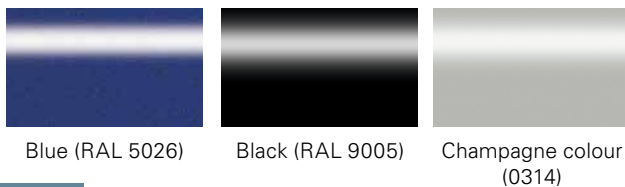

Blue (RAL 5026)    Black (RAL 9005)    Champagne colour (0314)

**Colour options**



Transparent red (0200)    White (RAL 9016)    Olive green (RAL 6003)

- 1/ Adapter to reduce seat depth by 130 mm available
- 2/ For small seat frame (490 mm) only
- 3/ Additional 50 mm of height adjustment
- 4/ Powered seat tilt
- 5/ Man. settable seat tilt
- 6/ according to ISO 7176-4:2008



More detailed information, technical data, references for transport or description of certain options can be found in the wheelchair's user manual.